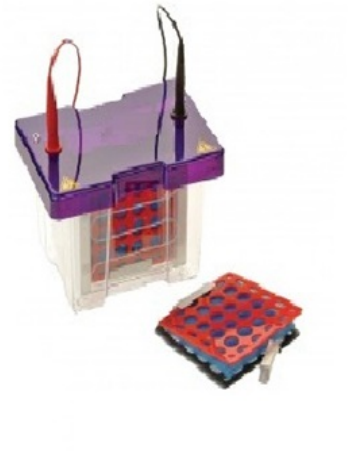




Vertical Sub Blot system AC

| Blotting

Details Info

Available in Mini, Mini Wide and Maxi formats, vertical Sub Blot systems are complete stand-alone units for submarine blotting. Each unit comprises an vertical tank and lid, electroblotting insert, and cassettes and fibre pads. The versatile design of the vertical tank and lid allows the unit to be easily adapted for vertical PAGE or IEF with capillary tube gels, using the relevant optional insert. Other features include:

- Multi-transfer capability for up to 3 Maxi gels, 6 Mini Wide and 12 Mini gels in the Maxi System; 3 Mini Wide and 6 Mini gels in the Mini Wide system; and 4 Mini gels in the Mini system
- Complete run-time flexibility from overnight to as short as a couple of hours
- Cool pack for environmentally friendly low-cost cooling during high-intensity transfers
- Open design, rigid cassettes maximise current transfer and eliminate 'shadow band' formation
- Cassettes colour-coded prevent polarity reversal to ensure transfer in the correct orientation

Ordering Information

ACSB10	OmniBlot Mini, 10 x 10cm Blotting System, including 4 cassettes
ACV10BI	OmniBlot Mini Insert - including 4 cassettes
ACV10BIEP	Wired Electrode Plate for ACV10BI (Replacement for Positive or Negative)
ACV10BI-HI	High Intensity Vertical Blot Mini insert - includes 2 cassettes and 8 fibre pads

SB10C	OmniBlot Mini Cassette
SB10F	Fibre pads - Pack of 8
SB10W	OmniBlot Mini Wide 20 x 10cm Blotting System, including 3 cassettes
ACV10WBI	OmniBlot Mini Wide Insert - including 3 cassettes
SB10WC	OmniBlot Mini Wide Cassette
SB10WF	Fibre pads - Pack of 6
SB20	OmniBlot Maxi, 20 x 20cm Blotting System, including 4 cassettes
ACV20BI	Sub Blot Maxi Insert - including 3 cassettes
SB20C	OmniBlot Maxi Cassette
SB20F	Fibre pads - Pack of 6
ACV10BI	omniBlot Mini Insert - including 3 cassettes
SB10C	omniBlot Mini Cassette
SB10F	Fibre pads - Pack of 8
SB10WC	Vertical Mini Wide, Cassette
ACV10WBI	VerticalSub Blot Mini Wide, Insert - including 3 cassettes
SB10WF	Fibre pads - Pack of 6
ACV20BI	Sub Blot Maxi Insert - including 3 cassettes
WAVEBI	Use Direct Replacment ACV20BI
SB20C	omni Blot Maxi Cassette
SB20F	Fibre pads - Pack of 6



Head Office :1600 AVENUE DE LORIMIER BUREAU, 384 MONTRÉAL (QUÉBEC) CANADA H2K3W8 USAOffice:128 W SUFFLOK AVE,CENTRAL ISLIP,11722 N Middle east & Africa Office : warehouse 9 /shed 9, Technology park P.O. Box 54806 RAK FZE , UAE



Email Us At:

sales@acculabusa.com middleeast@acculabusa.com



Phone Support: